

# Optimization of PCR Condition: The First Study of High Resolution Melting Technique for Screening of APOA1 Variance

*by* Udin Bahrudin

---

**Submission date:** 18-Sep-2019 08:57AM (UTC+0700)

**Submission ID:** 1174819902

**File name:** of\_PCR\_Condition\_The\_First\_Study\_of\_High\_Resolution\_Melting.pdf (1.45M)

**Word count:** 3933

**Character count:** 19391

## Optimization of PCR Condition: The First Study of High Resolution Melting Technique for Screening of *APOA1* Variance

Hesty Wahyuningsih,\* Ferdy K Cayami,\* Udin Bahrudin,\*† Mochamad A Sobirin,†‡ Farmaditya EP Mundhofir,\* Sultana MH Faradz\* and Ichiro Hisatome§

\*Center for Biomedical Research, Diponegoro University Faculty of Medicine, Semarang, Indonesia, †Department of Cardiology and Vascular Medicine, Diponegoro University Faculty of Medicine, Semarang, Indonesia, ‡Department of Pharmacology and Therapeutic, Diponegoro University Faculty of Medicine, Semarang, Indonesia, §Division of Regenerative Medicine and Therapeutics, Department of Genetic Medicine and Regenerative Therapeutics, Institute of Regenerative Medicine and Biofunction, Tottori University, Yonago 683-8503, Japan

### ABSTRACT

**Background** High resolution melting (HRM) is a post-PCR technique for variant screening and genotyping based on the different melting points of DNA fragments. The advantages of this technique are that it is fast, simple, and efficient and has a high output, particularly for screening of a large number of samples. *APOA1* encodes apolipoprotein A1 (apoA1) which is a major component of high density lipoprotein cholesterol (HDL-C). This study aimed to obtain an optimal quantitative polymerase chain reaction (qPCR)-HRM condition for screening of *APOA1* variance.

**Methods** Genomic DNA was isolated from a peripheral blood sample using the salting out method. *APOA1* was amplified using the RotorGeneQ 5Plex HRM. The PCR product was visualized with the HRM amplification curve and confirmed using gel electrophoresis. The melting profile was confirmed by looking at the melting curve.

**Results** Five sets of primers covering the translated region of *APOA1* exons were designed with expected PCR product size of 100–400 bps. The amplified segments of DNA were amplicons 2, 3, 4A, 4B, and 4C. Amplicons 2, 3 and 4B were optimized at an annealing temperature of 60 °C at 40 PCR cycles. Amplicon 4A was optimized at an annealing temperature of 62 °C at 45 PCR cycles. Amplicon 4C was optimized at an annealing temperature of 63 °C at 50 PCR cycles.

**Conclusion** In addition to the suitable procedures of DNA isolation and quantification, primer design and

an estimated PCR product size, the data of this study showed that appropriate annealing temperature and PCR cycles were important factors in optimization of HRM technique for variant screening in *APOA1*.

**Key words** *APOA1*; HRM; optimization

High resolution melting (HRM) is a post-PCR technique for variant screening and genotyping based on different melting points of DNA fragments.<sup>1</sup> The main principle of HRM is the integration of intercalating dye to double stranded DNA (dsDNA) and disintegration of the fluorescence dye when the dsDNA are separated to single stranded DNA (ssDNA). By using a real time measurement of the fluorescence level, a melting curve of fluorescence level against time is formed. The melting profile is mainly influenced by the length of DNA fragment, the GC content and heterozygosity.<sup>1–3</sup> Thus, the advantages of this technique are that it is fast, simple, efficient and high throughput result, particularly for screening of a large number of samples.<sup>1,3,4</sup>

*APOA1* is located in chromosome 11q23 and encodes apolipoprotein A1 (apoA1). *APOA1* spans about 2.2 kb and contains 4 exons.<sup>6</sup> ApoA1 composes up to 70% of HDL-C.<sup>7</sup> In the structure of HDL, apoA1 is a group of proteins arranged in an anti-parallel double-belt structure.<sup>8</sup> It is a cofactor for lecithin cholesterolacyl transferase (LCAT), substrate for scavenger receptor class B member 1 (SR-B1).<sup>9</sup> Deletions in the *APOA1* result in a decrease of HDL-C level to less than 40 mg/dL.<sup>10</sup> The monogenic disorder caused by *APOA1* mutation comprises about less than 5% of low HDL-C level cases,<sup>11</sup> but it accelerates the atherosclerosis process in contrast to other mutations in genes which encode the main protein of reverse cholesterol transport (RCT).<sup>12</sup>

This study aimed to obtain an optimal quantitative polymerase chain reaction (qPCR)-HRM condition for screening of *APOA1* variance. So far, it is the first study which screens variance in *APOA1* using HRM technique; therefore optimization of the technique is very

Corresponding author: Udin Bahrudin, MD, PhD  
bahrudin00@gmail.com

Received 2016 November 21

Accepted 2016 December 28

Abbreviations: *APOA1*, apolipoprotein class A member 1 gene; apoA1, apolipoprotein class A member 1 protein; dsDNA, double stranded DNA; HDL-C, high density lipoprotein cholesterol; HRM, high resolution melting; LCAT, lecithin cholesterolacyl transferase; qPCR, quantitative polymerase chain reaction; RCT, reverse cholesterol transport; SR-B1, scavenger receptor class B member1; ssDNA, single stranded DNA

important in order to get a reliable result.<sup>13</sup> In addition to a suitable procedure for DNA isolation and quantification, primer design and an estimated PCR product size, this study found that appropriate annealing temperature and PCR cycles were important factors in optimization of HRM technique for variant screening in *APOA1*.

## MATERIALS AND METHODS

### Primer design

The first step in this study was designing primers for *APOA1*. The referent sequences were taken from the database of UCSC Genome Bioinformatics Human Genome Browser Gateway (<http://genome.ucsc.edu/cgi-bin/hgGateway?redirect=manual&source=genome.ucsc.edu>) and Ensembl ([www.ensembl.org](http://www.ensembl.org)) with the accession numbers were GRCh38/hg38<sup>14</sup> and ENSG00000118137,<sup>6</sup> respectively. The primers were designed with PCR product sizes 100–400 bps using the *Lightcycler Primer Design Software* V.1.0.R.84 (Idaho Tech, Salt Lake City, UT). The concentration of primers both forward and reverse primers was 10 pmol/μL.

Based on the gene sequence database, the primers were designed from the start codon at exon 2. Exon 1 is an untranslated region, and thus it was not included in the analysis. The primers for amplicon 2 and amplicon 3 have a PCR product size of less than 250 bps as an ideal. They were located in intron with the purpose of covering the entire exon area. Exon 4 was divided into 3 amplicons, i.e. 4A, 4B and 4C in order to get optimum HRM results. The details are listed in Table 1. The

primers were designed to cover all exons from 2 to 4 and their flanking regions.

### DNA isolation and quantification

The genomic DNA was isolated using the salting out method. Then, the DNA concentration was quantified using the Nanovue Plus (GE Healthcare Bio Sciences AB, Buckinghamshire, UK). The DNA samples were diluted using distilled water to a concentration of 10 ng/μL. The genomic DNA from normal healthy donors was used in the optimization of qPCR-HRM. The criteria of normal healthy donors were that their HDL-C level > 40 mg/dL, cholesterol total level < 200 mg/dL, LDL-C level < 160 mg/dL, triglyceride level < 150 mg/dL, and their DNA sequence were similar to the reference sequence in the NCBI database. Informed consent for participation in this study was obtained. The investigation was approved by the Institutional Review Board of the Diponegoro University Faculty of Medicine and Dr. Ariadi Hospital (number 190/EC/FK/RSDK/2013) and conformed to the principles outlined in the Declaration of Helsinki.<sup>15</sup>

### HRM technique and its data analysis

*APOA1* was amplified using qPCR RotorGeneQ 5Plex HRM (Qiagen, Valencia, CA) with total volume 10 μL in each tube. Each tube consisted of 5 μL Type-It HRM Master Mix (Qiagen, Cat. No. 206542), 1 μL of forward and reverse primers each with a concentration of 10 pmol/μL, 2 μL of DNA with a concentration of 10 ng/μL and 2 μL of H<sub>2</sub>O. The curve analysis was done with

**Table 1. List of primers, nucleotide positions, PCR products size and their GC content**

Exon	Amplicon	Primer Sequences	Nucleotide Positions	Primer GC content (%)	Product size (bp)	Amplicon GC content (%)
2	2	F: GCCAGGCCCTTCTCTC	GRCh38/hg38: chr11: 116836271-116836287	64.7	116	59.0
		R: GGTTGGCTCCCTAGGTT	GRCh38/hg38: chr11: 116836370-116836386	58.8		
3	3	F: GGCCTGATCTGGGTCTC	GRCh38/hg38: chr11: 116836512-116836528	64.7	246	58.6
		R: TCATCCCACAGGCCTCT	GRCh38/hg38: chr11: 116836741-116836757	58.8		
4	4A	F: CCAGCCCTCAACCCCTTC	GRCh38/hg38: chr11: 116837258-116837274	64.7	290	60.3
		R: CCCTCTTGGAGCTCTGC	GRCh38/hg38: chr11: 116837531-116837547	58.8		
4	4B	F: GCAGGAGGAGATGGAGCTCTA	GRCh38/hg38: chr11: 116837485-116837505	57.1	327	69.7
		R: CCTTGGCGGAGGTCTCTC	GRCh38/hg38: chr11: 116837795-116837811	64.7		
4	4C	F: TCAGCGAGAAGGCCAAG	GRCh38/hg38: chr11: 116837769-116837785	58.8	273	63.4
		R: AGGAAGTCCCTGCTCCA	GRCh38/hg38: chr11: 116838025-116838041	58.8		

bp, base pair; F, forward primer; R, reverse primer; GC, guanine and cytosine.

Rotor Gene–Pure Detection version 2.1.0 (build 9) Software (Qiagen). The PCR product was visualized with the amplification curve and confirmed with 1.5% agarose gel for 30 minutes. The melting profile was analyzed through melting curve.

## RESULTS

Optimization of *APOA1* was done in exon 2 to 4 which were divided into 5 amplicons. They were amplified by using 5 pairs of primers. Because of its length, exon 4 was divided into 3 amplicons, *i.e.* 4A, 4B and 4C. Exon 1 were excluded. At 60 °C annealing temperature and 40 PCR cycles, the HRM result of amplicon 2, 3, 4B and 4C showed a crossing point less than cycle 30<sup>th</sup> and single peak of melting curve (Fig. 1A). Clear bands matched with the estimated PCR product sizes were seen in the agarose gel. However, double peaks of melting curve (Fig. 1B) and double bands on the agarose gel (Fig. 1C) were seen for amplicon 4A and 4C.

In order to get better HRM curves which represent correct product for amplicon 4A and 4C, we modified the annealing temperature and number of PCR cycle. For amplicon 4A, we tried to increase the temperature to 62 °C and 40 cycles. Although the PCR result showed one correct band in gel electrophoresis, the relative thin band and the amplification curve pattern (Fig. 2) indicated that PCR could be further optimized. PCR with temperature 62 °C and 45 cycles produced better results, *i.e.* the amplification curve increased after the 30<sup>th</sup> cycle (Fig. 3A), a single peak of melting curve (Fig. 3B), and a thick band matched with the estimated PCR product size was seen in the agarose gel (Fig. 3C). However, a smear was seen at the upper end of the estimated PCR product size. Further increasing the temperature to 63 °C and 50 cycles resulted in no amplification in the raw material of qPCR-HRM (data not shown).

Optimization of amplicon 4C was done with higher annealing temperature and PCR cycle than that in Fig. 2, *i.e.* 63 °C and 50 times, respectively. This condition resulted in an increase of amplification curve after the 30<sup>th</sup> cycle (Fig. 4A), a single peak of melting curve (Fig. 4B), and a clear single band matched with the estimated PCR product size in the agarose gel (Fig. 4C).

The summary of optimum qPCR-HRM condition of amplicons 2, 3, 4A, 4B and 4C are listed in Table 2.

## DISCUSSION

It has long been known that the important factors in the optimization of HRM technique for genetic variant screening are a suitable procedure of the DNA isolation and quantification, primer design and an estimated PCR product size,<sup>13</sup> and data from this study clearly supports



the idea that appropriate annealing temperature and PCR cycles were also important factors in optimization of HRM technique for variant screening in *APOA1*.

The primers covering exons 2, 3 and 4 were designed using a well-established software, the *Lightcycler Primer Design Software* V.1.0.R.84 (Idaho Tech). The maximum PCR product size was 327 bps in exon 4B which was qualified for maximum PCR product size to achieve 100% sensitivity and 100% specificity of the HRM technique.<sup>1, 13</sup> The primers of amplicons 2 and 3 were ideal primers with PCR product sizes of less than 250 bps,<sup>13</sup> while the size of amplicons 4A and 4C were 290 bps and 273 bps, respectively. These primer designs and the estimated PCR product size were clearly important in the optimization of PCR-HRM technique for variant screening in *APOA1*. In addition to target length, primer sequences should be carefully checked for known polymorphisms. Failed identification of any SNP in the primer sequence will reduce the success rate of HRM analysis since the HRM melting profile is influenced by DNA heterozygosity.<sup>1-3</sup>

Initially, the experiment used annealing temperature 60 °C and 40 PCR cycles and it showed that amplicons 2, 3, 4B and 4C have crossing points of less than 30 (Fig. 1A). It means that the sample was adequate and the amplification was efficient.<sup>13</sup> On the other hand, amplicon 4A had a crossing point more than 30 cycles, meaning that the sample may not be adequate and the amplification may not efficient.<sup>13</sup> Amplicons 4A and 4C had double bands on agarose gel. These could be caused by unspecific PCR products, and thus required an increase in annealing temperature to obtain optimal qPCR-HRM conditions. Other possible causes for results of amplicon 4A which had a crossing point more than 30, double peak melting curve and double band on agarose gel were less efficient DNA amplification, a low annealing temperature leading to unspecific PCR products, and a CG rich content (60.3%) in the DNA sequence.<sup>13</sup>



As shown in Fig. 2, an annealing temperature of 62 °C and 40 PCR cycles also failed to show an optimal condition for amplicon 4C. Thus, further increases of annealing temperature and PCR cycle, *i.e.* 63 °C and 50 cycles were performed to get a single peak of HRM curve and single band on agarose gel (Fig. 4).

Although a single peak of HRM curve could be obtained for amplicon 4A by increasing the annealing temperature from 60 to 62 °C, the amplification curve did not reach a plateau pattern (Fig. 2) suggesting that the optimum temperature was obtained but an additional PCR cycle was needed. Although a PCR with a temperature of 62 °C and 45 cycles produced better results, a thin smear remained at the upper the estimated PCR

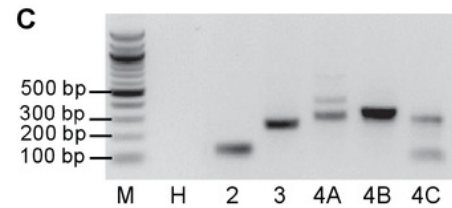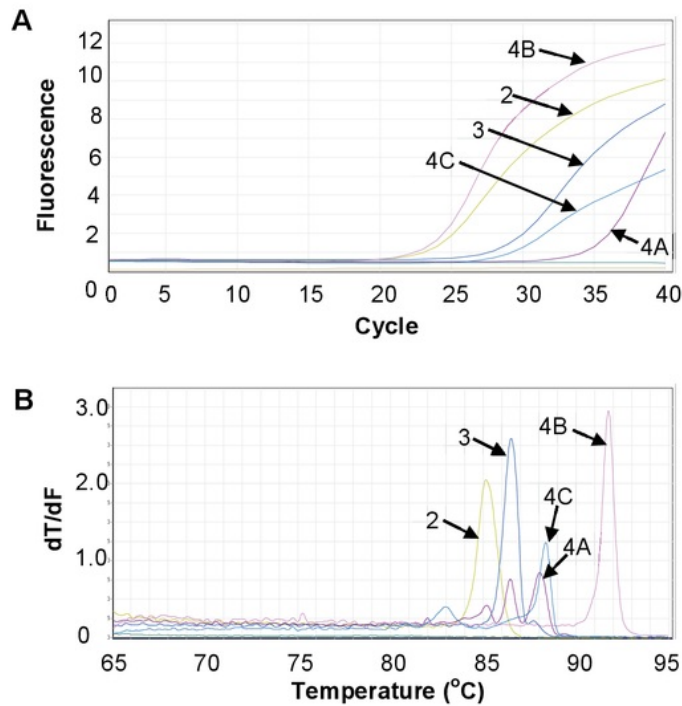


**Fig. 1.** The result of optimization of *APOA1* using annealing temperature 60 °C and 40 times of PCR cycles.

(A) qPCR amplification curve of amplicons 2, 3, 4B and 4C showed crossing point less than 30, while amplicon 4A showed that of more than 30. (B) HRM curve of amplicons 2, 3 and 4B had single peak, while amplicon 4A and 4C had double peak of melting. (C) Amplicons 2, 3 and 4B had clear single band on the agarose gel matched with the estimated PCR product size, but amplicons 4A and 4C had nonspecific band; the temperature of qPCR was set for starting at 65 °C. bp, base pair; dT/dF, delta temperature per delta fluorescence; H, H<sub>2</sub>O; HRM, high resolution melting; M, molecular weight marker; qPCR, quantitative polymerase chain reaction; 2, 3, 4A, 4B, and 4C are names of amplicons.

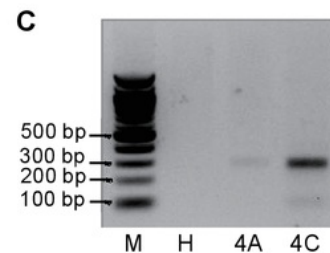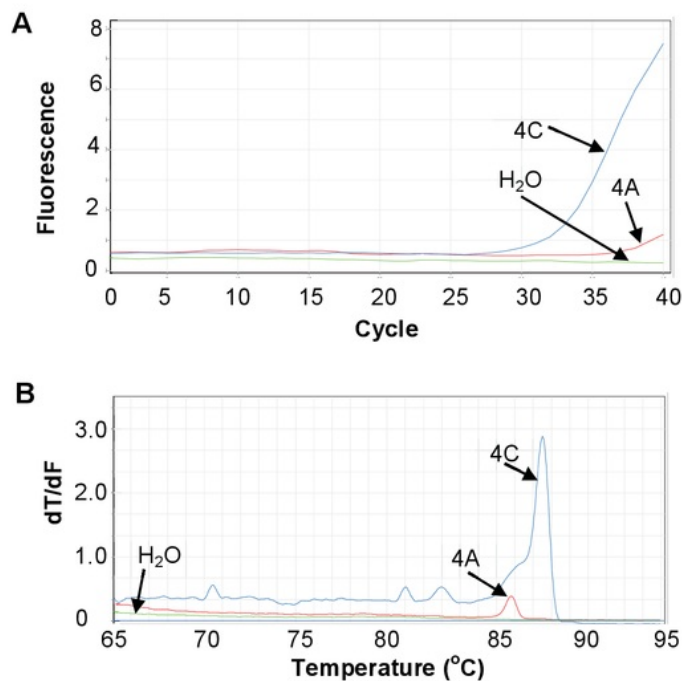


**Fig. 2.** The result of optimization of amplicon 4A and 4C *APOA1* using annealing temperature 62 °C and 40 times of PCR cycles.

(A) qPCR amplification curve of amplicons 4A and 4C showed crossing point more than 30. (B) HRM curve of amplicons 4A and 4C had single peak of melting. (C) Amplicon 4A had clear single band on the agarose gel matched with the estimated PCR product size, but amplicon 4C had a nonspecific band. bp, base pair; dT/dF, delta temperature per delta fluorescence; H, H<sub>2</sub>O; HRM, high resolution melting; M, molecular weight marker; qPCR, quantitative polymerase chain reaction; 4A and 4C are names of amplicons.

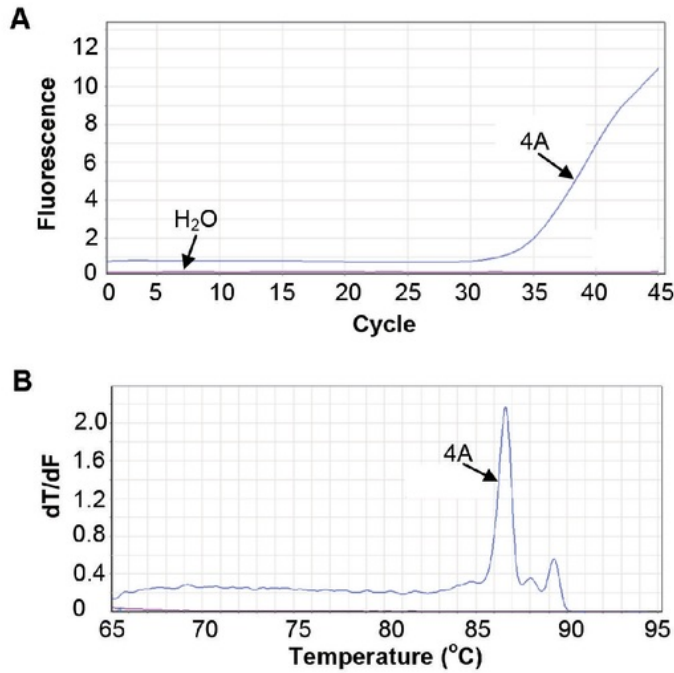

**Fig. 3.** The result of optimization of amplicon 4A *APOA1* using annealing temperature 62 °C and 45 times of PCR cycles.

(A) qPCR amplification curve of the amplicons showed that amplification was started at the 30<sup>th</sup> cycle. (B) Its HRM curve had single peak of melting. (C) The amplicon had a clear single band on the agarose gel matched with the estimated PCR product size, with very thin smear at the upper of the correct band. bp, base pair; dT/dF, delta temperature per delta fluorescence; H, H<sub>2</sub>O; HRM, high resolution melting; M, molecular weight marker; qPCR, quantitative polymerase chain reaction; 4A is a name of an amplicon.

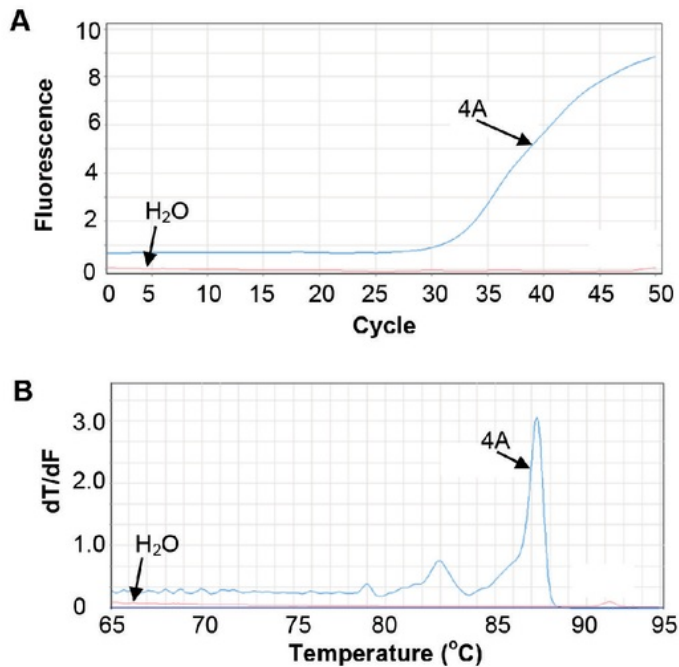

**Fig. 4.** The result of optimization of amplicon 4C *APOA1* using annealing temperature 63 °C and 50 times of PCR cycles.

(A) qPCR amplification curve of the amplicon 4C showed crossing point more than 30. (B) Its HRM curve had single peak of melting. (C) The amplicon had a clear single band on the agarose gel matched with the estimated PCR product size. bp, base pair; dT/dF, delta temperature per delta fluorescence; H, H<sub>2</sub>O; HRM, high resolution melting; M, molecular weight marker; qPCR, quantitative polymerase chain reaction; 4A is a name of an amplicon.

**Table 2. Profile of qPCR-HRM of *APOA1***

Amplicon	Program	Temperature (°C)	Time	Cycle
2, 3 and 4B	Initial denaturation	95	5 min	1
	Denaturation	95	10 s	40
	Annealing	60	10 s	
	Elongation	72	10 s	
	HRM	65–95 °C, 0.1 °C temperature increment Wait for 2 s for each steps afterwards		1
	HRM final step	40	1 min	1
	4A	Initial denaturation	95	5 min
Denaturation		95	10 s	45
Annealing		62	20 s	
Elongation		72	30 s	
HRM		65–95 °C, 0.1 °C temperature increment Wait for 2 s for each steps afterwards		1
HRM final step		40	1 min	1
4C		Initial denaturation	95	5 min
	Denaturation	95	10 s	50
	Annealing	63	10 s	
	Elongation	72	10 s	
	HRM	65–95 °C, 0.1 °C temperature increment Wait for 2 s for each steps afterwards		1
	HRM final step	40	1 min	1

HRM, high resolution melting; qPCR, quantitative polymerase chain reaction.

product size (Fig. 3), and further increasing of temperature to 63 °C and 50 cycles resulted in no amplification detected in the raw material of qPCR-HRM. Thus, at this point, considering that the amplification curve after the 30<sup>th</sup> cycle and a single peak of the melting curve could be achieved by using an annealing temperature of 62 °C and 45 cycles, this method could be accepted.

Our results indicated that there are primers which have similar PCR conditions for HRM, and thus facilitate multiple primers for a qPCR-HRM running. In contrast, some primers may need different conditions to achieve optimum conditions as in *APOA1*. To evaluate whether methods of *APOA1* variant screening showed in this study are appropriate for detection of mutant *APOA1*, further qPCR-HRM analysis using either genomic DNA or plasmid DNA harboring *APOA1* mutation is warranted.

In conclusion, in addition to a suitable procedure of DNA isolation and quantification, primer design and an estimated PCR product size, the data in this study showed that appropriate annealing temperature and PCR cycle were important factors in optimization of HRM technique for variant screening in *APOA1*.

**Acknowledgments:** This study was supported by a research grant from Diponegoro University (UB and MAS).

*The authors declare no conflict of interest.*

## REFERENCES

- 1 Reed GH, Kent JO, Wittwer CT. High-resolution DNA melting analysis for simple and efficient molecular diagnostics. *Pharmacogenomics*. 2007;8:597-608. PMID: 17559349.
- 2 AB Applied Biosystems [Internet]. CA: Life Technologies Corporation. A guide to high resolution melting (HRM) analysis. c2010-[cited 2016 Nov 19]. Available from: [http://www3.appliedbiosystems.com/cms/groups/mcb\\_marketing/documents/generaldocuments/cms\\_070283.pdf](http://www3.appliedbiosystems.com/cms/groups/mcb_marketing/documents/generaldocuments/cms_070283.pdf)
- 3 Wittwer CT. High-resolution DNA melting analysis: advancements and limitations. *Hum Mutat*. 2009 Jun;30:857-9. PMID: 19479960.
- 4 Vali Z, Raz A, Bokharaei H, Nabavi M, Bermanian MH, Yazdi MS, et al. Development of a High-resolution Melting Analysis Method Based on SYBR Green-I for rs7216389 Locus Genotyping in Asthmatic Child Patients. *Avicenna J Med Biotech*. 2014;6:72-80. PMID: 24834309.
- 5 Shanker J, Perumal G, Rao VS, Khadrinarasimhiah NB, John S, Hebbagodi S, et al. Genetic studies on the APOA1-C3-A5 gene cluster in Asian Indians with premature coronary artery disease. *Lipids Health Dis*. 2008;7:33. PMID: 18801202. PMCID: 2556320.
- 6 Ensembl [Internet]. Hinxton (UK): EMBL-EBI and WTSI. Gene Sequence Information for ENSG00000118137 [cited 2016 Nov 19] Available from: [http://asiaensemblorg/Homo\\_sapiens/Gene/Sequence?g=ENSG00000118137;r=11:116835751-116837950](http://asiaensemblorg/Homo_sapiens/Gene/Sequence?g=ENSG00000118137;r=11:116835751-116837950).
- 7 Rader DJ. Molecular regulation of HDL metabolism and function: implications for novel therapies. *J Clin Invest*. 2006 Dec;116:3090-100. PMID: 17143322. PMCID: 1679714.
- 8 Wasan KM, Brocks DR, Lee SD, Barrable KS, Thornton SJ. Impact of lipoproteins on the biological activity and disposi-

- tion of hydrophobic drugs: implications for drug discovery. *Nat Rev Drug Discov.* 2008;7:84-99. PMID: 18079757.
- 9 Wu Z, Wagner MA, Zheng L, Parks JS, Shy JM 3rd, Smith JD, et al. The refined structure of nascent HDL reveals a key functional domain for particle maturation and dysfunction. *Nat Struct Mol Biol.* 2007 Sep;14:861-8. PMID: 17676061.
- 10 Rye KA, Bursill CA, Lambert G, Tabet F, Barter PJ. The metabolism and anti-atherogenic properties of HDL. *J Lipid Res.* 2009;50 Suppl:S195-200. PMID: 19033213. PMCID: 2674714.
- 11 Kiss RS, Kavaslar N, Okuhira K, Freeman MW, Walter S, Milne RW, et al. Genetic etiology of isolated low HDL syndrome: incidence and heterogeneity of efflux defects. *Arterioscler Thromb Vasc Biol.* 2007 ;27:1139-45. PMID: 17303779.
- 12 Hovingh GK, de Groot E, van der Steeg W, Boekholdt SM, Hutten BA, Kuivenhoven JA, et al. Inherited disorders of HDL metabolism and atherosclerosis. *Curr Opin Lipidol.* 2005;16:139-45.
- 13 Lightcycler@480 Real Time PCR System [Internet]. Mannheim (Germany): Roche Diagnostics. Technical note no.1, high resolution melting: optimization strategies. c2008- [cited 2016 Nov 19]. Available from: [https://lifescience.roche.com/wcsstore/RASCatalogAssetStore/Articles/05382238001\\_11.08.pdf](https://lifescience.roche.com/wcsstore/RASCatalogAssetStore/Articles/05382238001_11.08.pdf).
- 14 UCSC Genome Browser [Internet]. CA: UC Santa Cruz. Gene Sequence Information for NM\_000039 [cited 2016 Nov 19]. Available from: [https://genome.ucsc.edu/cgi-bin/hgTracks?db=hg38&lastVirtModeType=default&lastVirtModeExtraState=&virtModeType=default&virtMode=0&nonVirtPosition=&position=chr11%3A116835751-116837950&hg-sid=568378555\\_QBOqVJC0HYbY6GNCEKvUaJSsm0dZ](https://genome.ucsc.edu/cgi-bin/hgTracks?db=hg38&lastVirtModeType=default&lastVirtModeExtraState=&virtModeType=default&virtMode=0&nonVirtPosition=&position=chr11%3A116835751-116837950&hg-sid=568378555_QBOqVJC0HYbY6GNCEKvUaJSsm0dZ).
- 15 World Medical Association Declaration of Helsinki. Recommendations guiding physicians in biomedical research involving human subjects. *Cardiovasc Res.* 1997;35:2-3.

# Optimization of PCR Condition: The First Study of High Resolution Melting Technique for Screening of APOA1 Variance

## ORIGINALITY REPORT

1 %	2 %	2 %	%
SIMILARITY INDEX	INTERNET SOURCES	PUBLICATIONS	STUDENT PAPERS

## PRIMARY SOURCES

1	Bahrudin, U.. "Ubiquitin-Proteasome System Impairment Caused by a Missense Cardiac Myosin-binding Protein C Mutation and Associated with Cardiac Dysfunction in Hypertrophic Cardiomyopathy", Journal of Molecular Biology, 20081226 Publication	1 %
2	<a href="http://www.jstage.jst.go.jp">www.jstage.jst.go.jp</a> Internet Source	1 %

Exclude quotes	On	Exclude matches	< 1 %
Exclude bibliography	On		